





IOL Removal Set

1-705 Osher Angled IOL Cutter

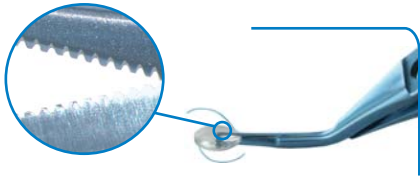

Micro serrations grip the IOL as the scissors cut. IOL is kept stable and does not slip in the blades of the scissors.

Suitable for Silicone and Acrylic IOL's

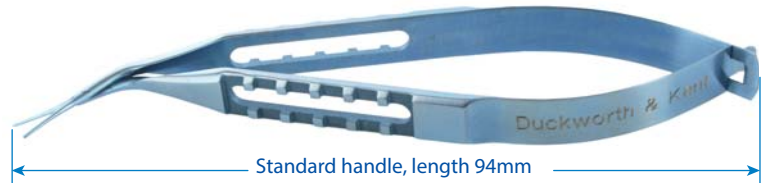


Standard handle, length 94mm

- Blunt tips with micro serrated blades
- Cut length 8mm
- Tip to pivot length 9.5mm

video available at
www.duckworth-and-kent.com/videos

2-796N DK IOL Removal Forceps

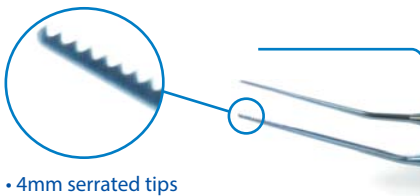

- 4mm serrated tips
- 30° angled shafts, tip to angle length 14mm



Standard handle, length 82.5mm


1-805 Osher Haptic Cutter



- 20 gauge, 20mm long shaft
- Squeeze action activates jaws to sever haptic



Round squeeze handle, length 145mm

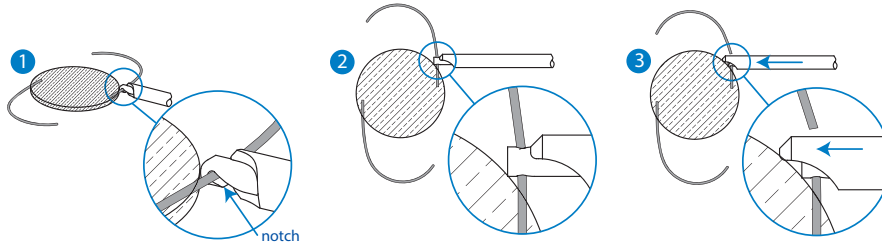

Tips of the the Osher Haptic Cutter have a notch that fits over and grasps the haptic (fig1). When the handle of the Osher Haptic Cutter is squeezed (fig 2 - 3) the outside tube of the tip moves forwards and severs the haptic.

www.duckworth-and-kent.com

7 Marquis Business Centre
Royston Road, Baldock
Herts SG7 6XL England

Tel: +44 (0)1462 893254
Fax: +44 (0)1462 896288
Email: info@duckworth-and-kent.co.uk



© March 2009 Duckworth & Kent
Revised 01.07.2016

D&K® is a registered trademark. All other brand names are trademarks or registered trademarks of their respective owners. All schematic line drawings, photographs and copy in this leaflet are fully protected by copyright. No part of this leaflet may be reproduced in any form without prior written permission. We reserve the right to make changes at any time, without notice, in product specifications and availability. Descriptive, typographic, or photographic errors are subject to correction. Name(s) of instruments are often comprised of surgeon's name, combination of surgeons' names or by the category of the instrument.

at the Leading Edge