




Suggested Khaw Glaucoma Set



A set of instruments amalgamating those designed by Professor Khaw to assist the surgeon in performing trabeculectomy using the "safe surgery system".

9-576 Khaw Standard Glaucoma Surgery Speculum

- 14.5mm solid blades
- Angled to rest temporarily
- Adjustable with thumb screw

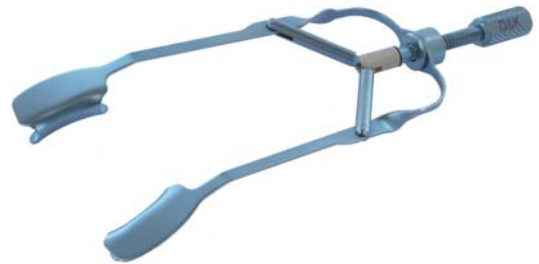


8-616 Ogawa Infusion Cannula

- 20 gauge, thin wall, 0.9mm diameter x 4mm length
- Straight shaft, 45° bevelled tip
- Grooves in cannula body ensure non-traumatic securing of cannula
- Silicone tubing and luer fitting supplied




2-686 Khaw Small Conjunctival Clamp

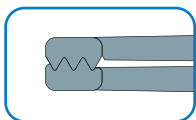

- Holds conjunctiva securely
- Holds and protects the conjunctival edge during antimetabolite application
- Particularly during fornix based conjunctival incisions
- Single handed action
- Tip width 4mm
- Overall length 74.5mm




2-687 Khaw Large Conjunctival Clamp

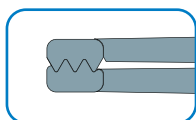

- Holds conjunctiva securely
- Holds and protects the conjunctival edge during antimetabolite application
- Particularly during fornix based conjunctival incisions
- Single handed action
- Tip width 12mm
- Overall length 74.5mm




7-101 Khaw Small Descemet Membrane Punch

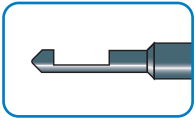


- Designed to punch 0.5mm x 0.3mm
- Particularly suited when small sclerostomy required (e.g. thin sclera, small scleral flap)
- Punch action can be repeated to create larger sclerostomy
- Can be used with a short scleral tunnel incision
- Squeeze action handle activates shaft to punch

7-102 Khaw Descemet Membrane Punch

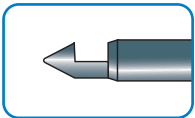

- Designed to punch 0.75mm x 0.5mm
- Can be used with a short scleral tunnel incision
- Squeeze action handle activates shaft to punch
- Round squeeze handle, length 131mm


2-502N Khaw Transconjunctival Adjustable Suture Control Forceps

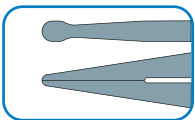

- 5.0mm highly polished tying platforms with a flared tip
- For massaging and adjusting intraocular pressure to desired level
- Adjustable Suture Control technique
- Straight shafts
- Standard handle, length 85mm

DK7003 Sterilising Tray

- External dimensions: 17.5cm x 11cm x 2.6cm
- Internal dimensions: 16.5cm x 10cm x 1.8cm (without silicone mat)
- One silicone mat

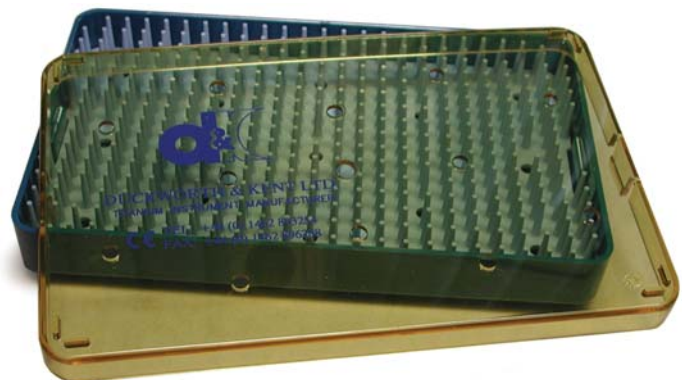

www.duckworth-and-kent.com

7 Marquis Business Centre
Royston Road, Baldock
Herts SG7 6XL England

Tel: +44 (0)1462 893254
Fax: +44 (0)1462 896288
Email: info@duckworth-and-kent.co.uk



© March 2011 Duckworth & Kent
12.06.16

DK&K is a registered trademark. All other brand names are trademarks or registered trademarks of their respective owners. All schematic line drawings, photographs and copy in this leaflet are fully protected by copyright. No part of this leaflet may be reproduced in any form without prior written permission. We reserve the right to make changes at any time, without notice, in product specifications and availability. Descriptive, typographic, or photographic errors are subject to correction. Name(s) of instruments are often comprised of surgeon's name, combination of surgeons' names or by the category of the instrument.

at the Leading Edge